[PROF]

RESPIRATOR WITH VALVE

60156490

Instruction manual (Original instructions)	GB	1
Käyttöohje (Alkuperäisten ohjeiden käännös)	FI	2
Bruksanvisning (Översättning av originalinstruktionerna)	SE	3
Bruksanvisning (Oversettelse av de opprinnelige instruksjonene)	NO	4
Kasutusjuhend (Tõlgitud originaal juhendist)	EE	5
Instrukcija (Tulkojums no oriģinālvalodas)	LV	6
Naudojimo instrukcijos (Vertimas originali instrukcija)	LT	7
Instrukcja obsługi (Oryginalna instrukcja)	PL	8

USE INSTRUCTIONS


1. Failure to follow all instructions and limitations can seriously Reduce the effectiveness of this respirator and could lead to Illness, injury or death.
2. A properly selected respirator is essential. Before occupational use, the wearer must be trained by the employer in the correct use of the respirator in accordance with applicable safety and health standards.
3. This product does not supply oxygen. Use only in adequately ventilated areas containing sufficient oxygen to support life.
4. Discard the respirator and replace with a new one if:
 - a. excessive clogging of the respirator causes breathing difficulty.
 - b. the respirator becomes damaged.
5. Leave the contaminated area if dizziness, irritation or other Distress occurs.
6. Keep respirators in the display box away from direct sunlight until use.
7. Before wearing the mask respirator users must check whether there is crack or spot around the welding spot. if yes, the mask must be discarded at once and can not be used again.
8. Keep the unused respirators in their closed box and store in a dry non-contaminated area between -30 °C and +50 °C. Storage life is 3 years from date of manufacture.

USE LIMITATIONS


1. Do NOT use the respirator or enter or stay in a contaminated area under the following circumstances:
 - a. atmosphere contains less than 19.5% oxygen.
 - b. if you smell or taste contaminant.
 - c. for protection against gases or vapors.
 - d. contaminants or their concentrations are unknown or immediately dangerous to life or health.
 - e. concentrations of contaminants exceed maximum use concentrations in applicable OSHA standards or applicable government regulations or 12 (FFP2) times the PEL (Permissible Exposure Limit), whichever is lower.
 - f. For sandblasting, paint-spray operations and asbestos.
 - g. Use in explosive atmospheres.
2. Do NOT modify or misuse the respirator.
3. Do NOT use the respirator with facial hair or any other conditions that may prevent a good face-seal.
4. NR means non reuse. Maximal application time or maximal particle exposure one single shift.
5. D-dolomite, The particle filtering half masks have passed the storage test with dolomite.
6. FFP means Filtering Facepiece.
7. No filtering respiratory protection in containers, tanks, channels.

FITTING INSTRUCTIONS

1. Hold the respirator in hand with the contoured section of the respirator pointing up. Allow headbands to hang freely.
2. Position the respirator under the chin, covering mouth and nose.
3. Pull the lower headband over the head, position below ears. Pull the top headband on the back of head above ears for best fit. Bend the nose clip around fingers to form nose shape.
4. Using both hands, form the nose clip around the nose and shift the headstraps to ensure a positive face seal.
5. Seal check:
 - To test the fit of a respirator without an exhalation valve, cup both hands over the respirator and exhale sharply.
 - To test the fit of a respirator with an exhalation valve, cup both hands over the mask and inhale sharply. A negative pressure should be detected inside the respirator.
 - If air flow is felt in the nose area, re-adjust/tighten the nose clip.
 - If air flow is felt around the edges of the respirator, re-position the respirator/headstrap to achieve a better fit.
6. Change the mask immediately if breathing becomes difficult or mask becomes damaged or distorted.
7. Change the respirator if a proper face seal can not be achieved.
8. Careful observance of these instructions is an important step in safe respirator use.

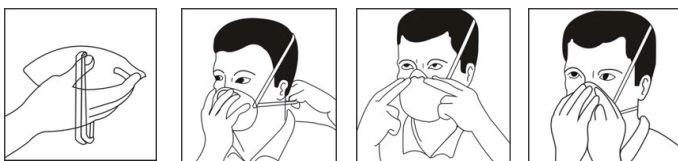See information
supplied by the
manufacturer

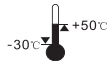Temperature
range of storage
conditions

Maximum relative
humidity of storage
conditions

WARNING

Conventional carton packing. Avoid damage to the package during transportation, and avoid pressed by heavy object, avoid moisture or sun exposure.

This respirator helps protect against certain particulate contaminants but does not eliminate exposure to or the risk of contracting disease or infection. Misuse may result in sickness or death. Hair, beards, sideburns and anomalies of the head under seal line will affect the seal of the respirator.

These filtering respiratory are in conformance with the provision of PPE regulation 2016/425 and, where such is the case, with the harmonized standard EN149:2001+A1:2009 is identical to the PPE which is the subject of EC type.

For all these issues if they impact the safety with the end users, please contact with the manufacture and BSI.

Declaration of conformity: <https://www.byggmakker.no/produkt/STOVMASKE-FFP2-MED-VENTIL-PROF/7070483142958>

**BSI Group - NB 2797
The Netherlands BV, Say Building, John M Keynesplein 9,
1066 EP, Amsterdam, Netherlands.**

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Työpajankatu 12, FI- 00580 Helsinki © Kesko 2024. Made in China.
Gamintojas • Producent • Kesko Corporation Building and technical trade,